
Rees Scientific Corporation – Centron Environmental Alarm Instructions for Telephoned Alarms

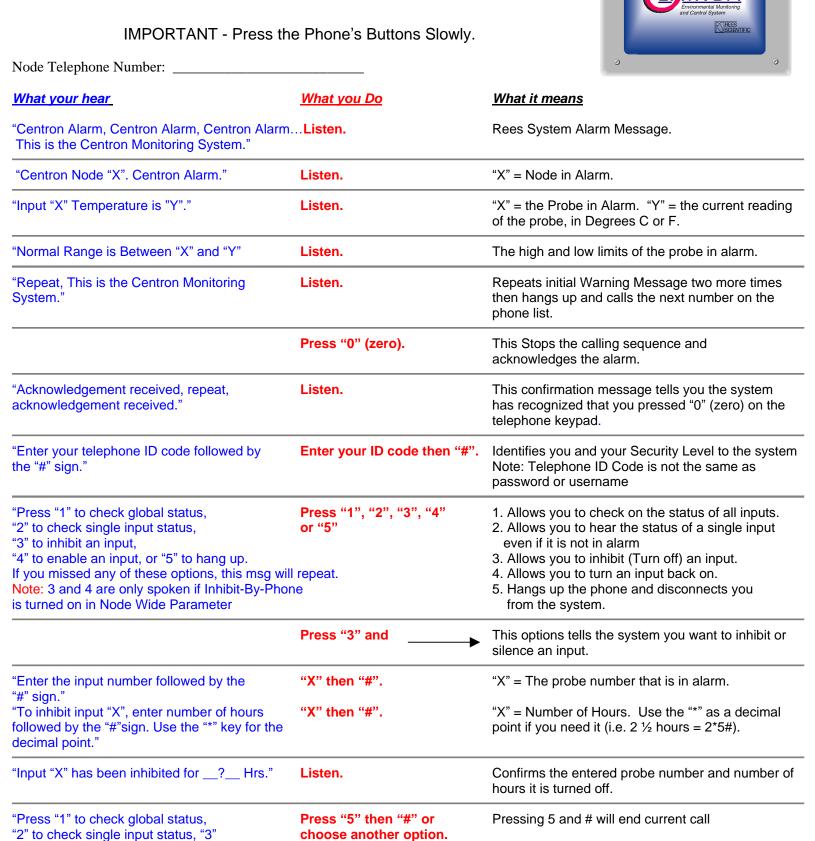

The Message will repeat itself. If more than one probe is in alarm, you would take care of them one at a time.

"5" to hang up."

to inhibit an input, "4" to enable an input,

Rees Scientific Corporation - Centron Environmental Alarm

BASIC TIPS

- Don't give out your pass code to anyone. With your code other people are able to make changes to the system (at your
- security level). That will make it look like you did it.
- 3. Utilizing telephone option number 1 or 2, you can call back into the system at a later time and check on the real-time status of the probe you inhibited (or any probe for that matter). This way you can see if it is still moving further out of range or coming back into range. (Use the Nodes telephone number on front page)
- 4. Always report any telephone activity with the Centron system to your supervisor. Documentation needs to be maintained on any interaction with the system. (i.e. the why of what you did.) The system records the who, what, and when of any interaction.
- 5. Notify your supervisor of any changes to your telephone number.

Press the Telephone's buttons SLOWLY.

6.	Know which probe number refers to which piece of equipment number below.	. Write the equipment you are responsible for with the probe

1.